

Fast analysis of powders and defects in additively manufactured components

Metal additive manufacturing has emerged as a promising new technology primarily for aerospace and biomedical component applications. Technologies such as powder-bed or directed energy deposition require the use of metal powders. The powder chemistry, size distribution, and sphericity play key roles for the beam-induced melting process. While established technologies exist to examine powder size distributions (for example, laser-diffraction-based techniques), in some instances, powders are now used with sub-micrometer sizes that are very difficult to assess with these techniques. Existing technologies to analyze powders cannot distinguish between different elements or chemistries.

Challenge

Despite the novel and advanced nature of additive manufacturing, many of the issues that have long troubled traditional manufacturing techniques, such as defects and porosity, also emerge in additively manufactured parts. Therefore, much focus is currently placed on determining the nature of defects in additively manufactured parts and their link to the processing parameters. Here again, a manual defect-by-defect analysis using electron microscopy cannot yield defect size and shape distributions with the sampling sizes that are required by many manufacturers.


Method

In this application note, we show how the automated analysis mode for the Thermo Scientific™ Explorer™ 4 Analyzer can be used to reveal powder particle size distributions and defect distributions. Not only can size and shape distributions be obtained, but the EDX detector can also be used to differentiate between powder particles or microstructural features with different chemistries. For this example, commercial Ti-6Al-4V powders were obtained from two vendors, hereafter referred to as Supplier A and Supplier B.


The Explorer 4 Analyzer enables statistical analyses of microstructures and particles with both morphology and chemistry filters. For example, Figure 1a shows a blend of pure titanium and Ti-6Al-4V powder particles. Both powder particle types are nearly perfectly spherical and their Z contrast is similar. The two powder particle types are therefore not distinguishable with backscattered electron (BSE) imaging. However, with the Explorer 4 Analyzer, a strategy can be devised to analyze either one of the two powder types.

Analysis


In the first step of the analysis, the system software separates all particles from the background using an image contrast threshold criterion. Once the scanning electron beam detects features with

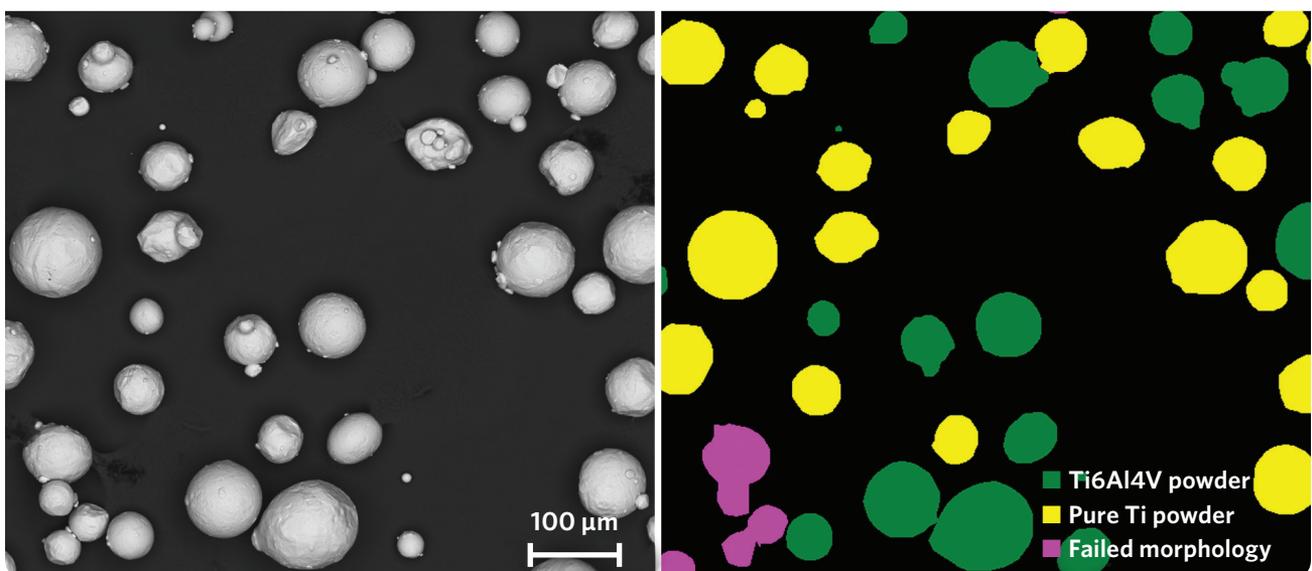


↑ **Figure 1.** Example of chemistry-specific particle analysis for Ti-6Al-4V powder using the Explorer 4 Analyzer.

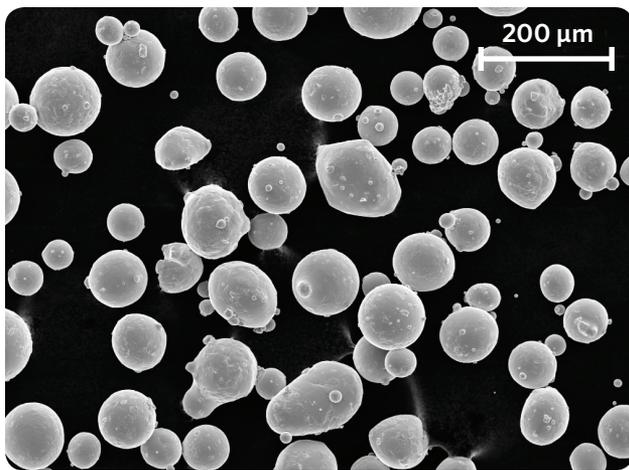


(a)

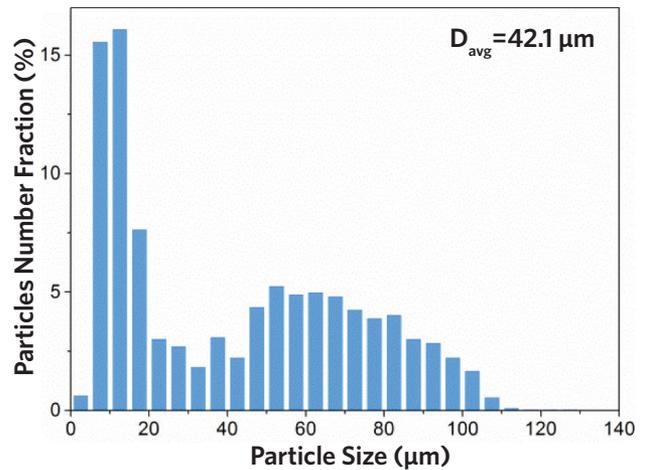

(b)

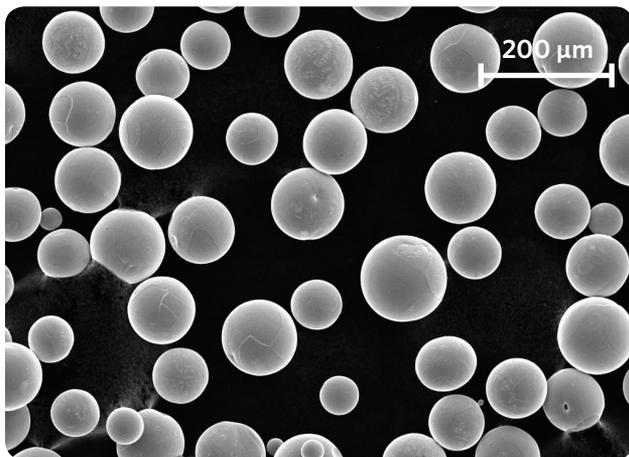

(c)

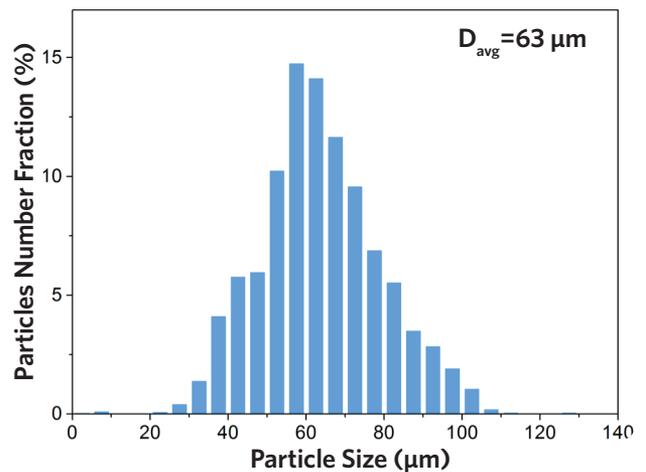


(d)

↑ **Figure 2.** Secondary electron SEM images and the corresponding particle size distributions measured by the Explorer 4 Analyzer for Supplier A (a, b) and Supplier B (c,d) Ti-6Al-4V powder.

image pixel intensities above the threshold, the instrument performs a short-duration (3 seconds in this work) EDX measurement of these features as a rough chemistry assessment. The rough chemistry information is then used for a rule set with morphology and chemistry filters.

In this example, the result of applying the filters is shown in Figure 1b. The yellow color-coded pure titanium powder was rejected by the chemistry composition filters, and the scanning electron beam skipped these features. Precise size and EDX chemistry analysis was then carried out for only the particles that passed chemistry and morphology filters, shown in Figure 1b as green and purple color-coded Ti-6Al-4V powder particles. The particle aspect ratio can be used as an additional filter. A threshold was set for the aspect ratio so that green color-coded spherical or nearly-spherical particles were separated from the purple color-coded irregular-shaped particles. This approach avoids complications from overlapping or otherwise agglomerated particles contributing to the particle size analysis. As a result of the chemistry and morphology rule set, only spherical Ti-6Al-4V powder particles were analyzed in detail.

Examples of automatically measured powder particle size distributions are shown in Figures 2 and 3. For the two size distributions shown in Figure 2, the software analyzed approximately 4,000 particles for each sample, corresponding to an analysis time of 20 minutes per sample. The two size distributions differ in the number of particles with sizes less than about 30 μm: the powder from Supplier A contains a large number of sub-30 μm particles, while the powder from Supplier B reveals a cut-off particle size of approximately 25 μm. The measurement time required to obtain the distributions can be reduced by limiting the particle count. Figure 3 shows the effects of particle count on the size distribution for counts of 1,000–10,000 particles from the same sample. The distributions show very little change for counts of more than 2,000 particles. The median particle size, in fact, only changed by approximately 1 micrometer upon extending the measurement time from 10 minutes to 50 minutes (i.e. from 2,000 particles to 10,000 particles). These measurements highlight the potential of the Explorer 4 Analyzer to analyze powder particles for sizedistribution and by chemistry with measurement times on the order of 10 minutes.

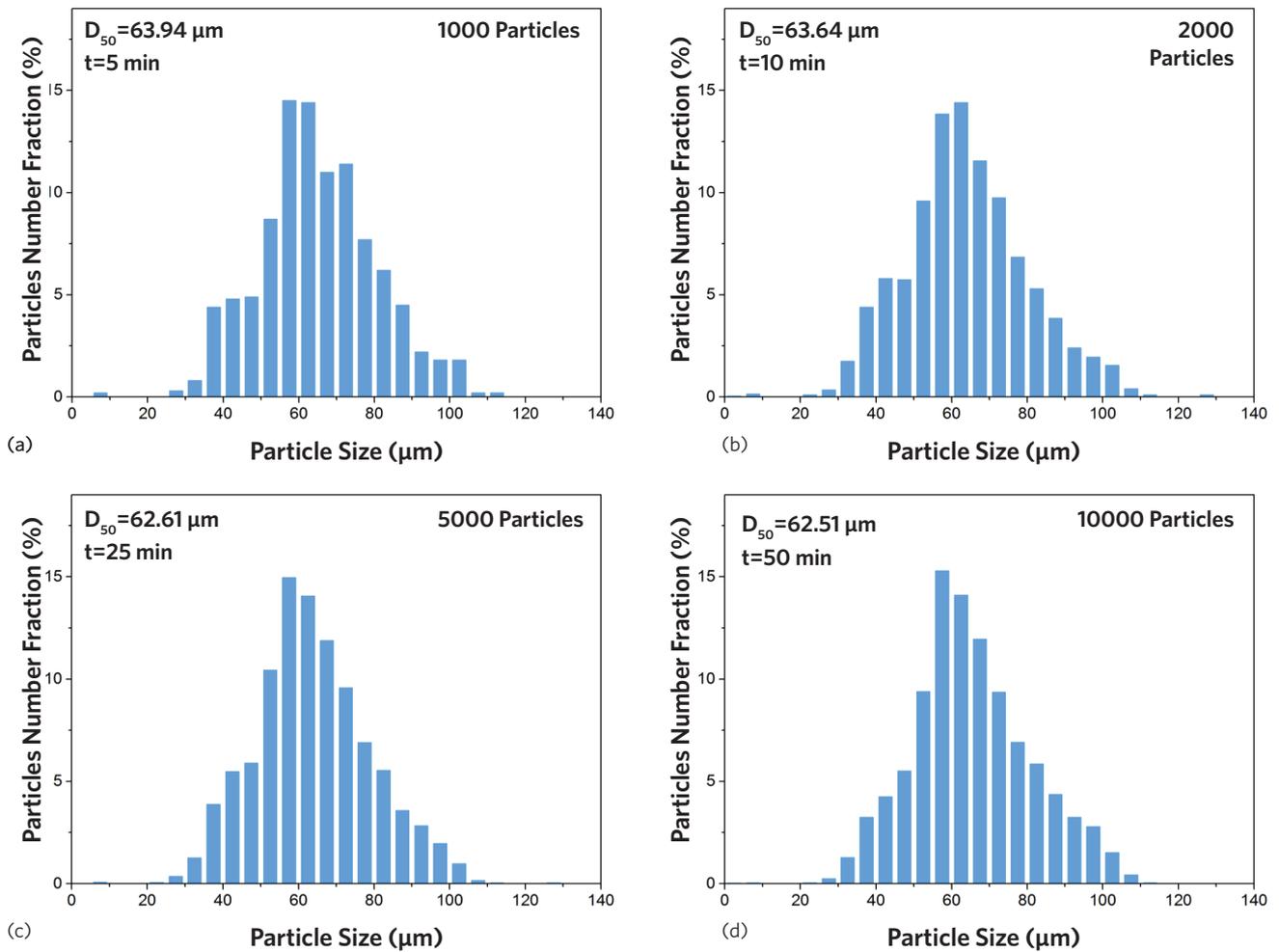


↑ **Figure 3.** Particle size distributions measured from Supplier B Ti-6Al-4V powder in the Explorer 4 Analyzer using: 1,000 particles (a), 2,000 particles (b), 5,000 particles (c), and 10,000 particles (d).

Defects


After parts are made with additive manufacturing, the focus of the analysis shifts from the powder pedigree to the microstructures and defects in the additively manufactured parts. Some defects, such as unmelted regions, are typically on the order of micrometers to hundreds of micrometers and can thus be identified even with light microscopy. Smaller defects, such as pores or inclusions, are detectable only with electron microscopy. The Explorer 4 Analyzer and its automated measurement capabilities are well suited to generate defect maps and size distributions over regions with dimensions of centimeters. Figure 4a is a photograph of an additively manufactured test sample. Figure 4b is an image from a cross-section through the test sample showing an unmelted region next to a spherical inclusion. Within the elongated, unmelted region, individual round powder particles can still be discerned. The cause of these defects is certainly of interest, but equally relevant is their distribution within the sample and their size distribution. Since the spherical inclusion differs in chemistry from the matrix and the unmelted region, the chemistry-based filtering mode of the Explorer 4 Analyzer allows for an analysis of size distributions only of inclusions.

An example of the system's capability to sample macroscopic regions for defects is shown in Figure 5. This is an analysis of large-scale cross-sectional areas for the additively manufactured test sample shown in Figure 4a. The 1 mm x 1 mm square area (colored pink in Figure 5a), resulted in pore and inclusion maps that are shown in Figure 5b as individual and combined maps. The software allows macroscopic regions to be defined for analysis, and the instrument then divides these regions into smaller frames for measurement purposes.

For example, an area of the main portion of the additively manufactured sample, shown with a black rectangle in Figure 5a, was analyzed for pore and inclusion volume fractions. The same analysis was applied to the thin section on top of the sample. The results of these measurements are shown in Figure 5c. For this particular sample, the pore volume fraction (taken as being equal to the area fraction, based on standard quantitative stereology methods) is 2.3% for the thick section and 0.25% for the thin section, while the fraction of inclusions is 8.6% for the thick section and 3% for the thin section.

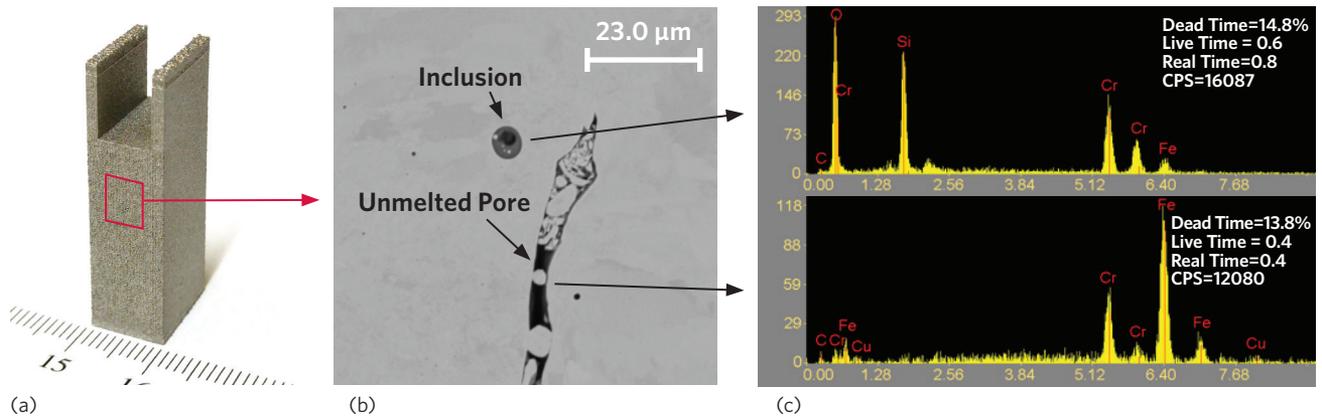

↑ **Figure 4.** (a) Photograph of the AM as-built test part, (b) Backscattered electron image of a typical unmelted pore and an inclusion in the AM as-built test part, (c) Examples of EDS spectra collected on a silicon oxide inclusion (top) and unmelted pore (bottom) during automatic feature analysis.

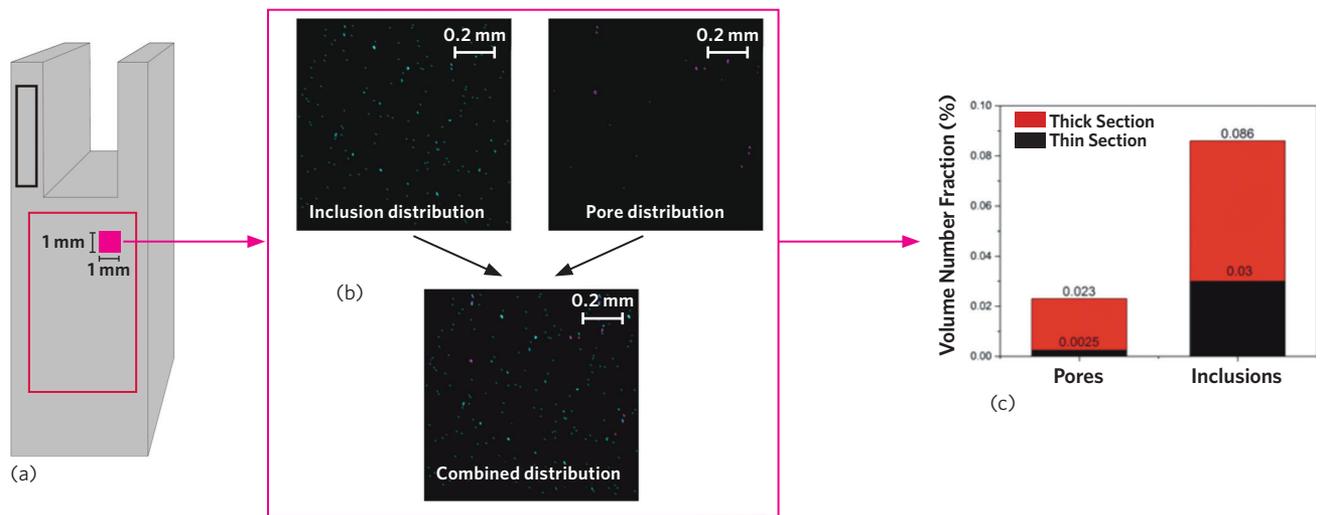

↑ **Figure 5.** (a) Schematic illustration of an additively manufactured test sample with analysis regions (pink box, red and black rectangles). (b) Measurement results for inclusion and pore distribution and combined image for the 1 mm x 1 mm area highlighted in pink in (a). (c) Measured overall volume fractions of pores and inclusions in the thick and thin sections of the AM as-built part (shown with red and black rectangles, respectively, in (a)).

Results

The Explorer 4 Analyzer allows for quick—approximately 10 minutes—analysis of powder particle size distributions and defect size distributions. Additionally, the system can analyze and differentiate based on chemical composition and morphologies. For additive manufacturing, key applications include the analysis of defects in microstructures of additively manufactured parts and powders.

Acknowledgments

This work was conducted under the FEI and University of Connecticut collaboration agreement at the Center for Advanced Microscopy and Materials Analysis with samples prepared at the UConn Additive Manufacturing Innovation Center.

World Headquarters
Phone +1 503 726 7500

FEI Europe
Phone +31 40 23 56000

FEI Japan
Phone +81 3 3740 0980

FEI Asia
Phone +65 6272 0050

FEI Australia
Phone +61 2 6173 6200

Learn more at FEI.com
ContactUs@FEI.com

For current certifications, please visit FEI.com/certifications

©2017. We are constantly improving the performance of our products—all specifications are subject to change without notice. FEI, the FEI logo, and Thermo Scientific Explorer 4 Analyzer are trademarks of FEI Company or its affiliates. All other trademarks belong to their respective owners.
AN0073-EN-04-2017